
C N I C P P T 0 N L Ι A L R E A R A I

To assure clinical success, the following minimum material thickness is recommended.

Preparation:


Anterior Crown

- 1.5 2.0 mm incisal reduction
- 1.0 1.5 mm labial and lingual reduction
- Round the internal line angles
- Chamfer margin

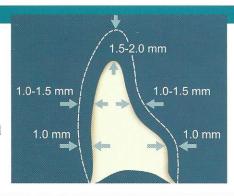


Posterior Crown


- 1.5 2.0 mm occlusal/incisal reduction
- 1.0 2.0 mm axial reduction
- Round the internal line angles
- Chamfer margin

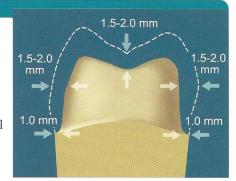

Preparation:

The optimal preparation is a shoulder or chamfered preparation with a circumferential step or chamfer which must be applied at an angle of >5° (horizontal). The angle of the preparation (vertical) should be 4° or larger. The inside angle of the shoulder preparation must be given a rounded contour.

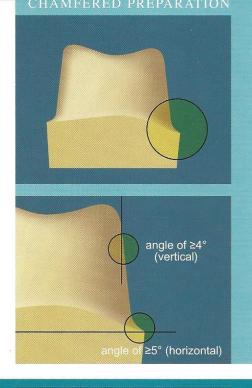

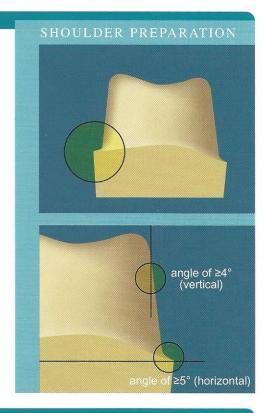

Contraindication:

We do not recommend a feather edge, as it does not provide adequate reduction for the porcelain buildup, or the trough (or gutter) shoulder because the outer most edge may not be detected when scanned.

